

Premium quality. Complete endo essentials.

NIC is the international Brand of Superline Technology.

Founded in 2002, Superline specializes in endodontic and orthodontic products. Its 10,000 m² facility combines advanced CNC manufacturing, precision technology and in-house R&D.; With over 20 patents, NIC serves dental professionals in more than 70 countries.



H-FILES



K-FILES



U-Blue File System



ISO PAPER POINTS



U-Series Paper Points • F1-F3



ISO GUTTA-PERCHA



U-Series Gutta-Percha • F1-F3

HAND INSTRUMENTATION



[Click to Visit Website](#)




H-Files & K-Files

Stainless steel hand files for canal preparation.

Stainless Steel


21 mm | 25 mm

H-Files

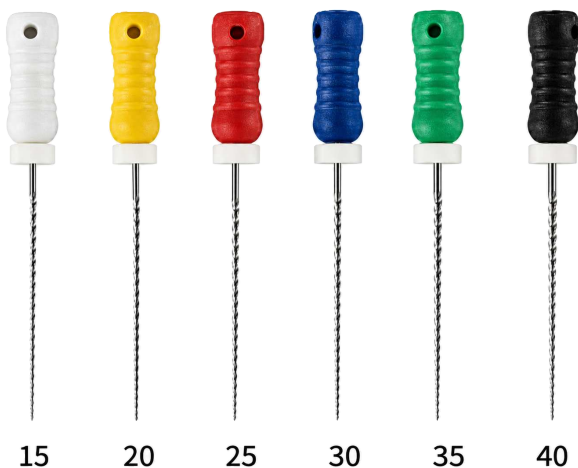

Designed for efficient debris removal and root canal enlargement. Their stainless steel construction and sharp cutting edges support controlled manual preparation, making them a practical choice for everyday endodontic treatment.

K-Files

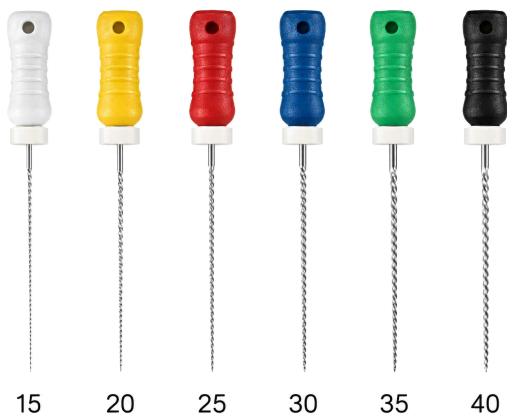

Designed for root canal exploration and manual preparation. They support controlled instrumentation & colour-coded handles simplify size identification during treatment.

ROTARY SHAPING



[Click to Visit Website](#)



U-Blue rotary files

MaxTech 6.0 NiTi alloy for sequential canal shaping.

Flexible NiTi technology.

Variable taper design.

Continuous clockwise rotation.




AVAILABLE LENGTHS

SX: 19 mm | S1-F5: 21 mm and 25 mm

NiTi alloy

360° rotation

Convex triangular cross-section

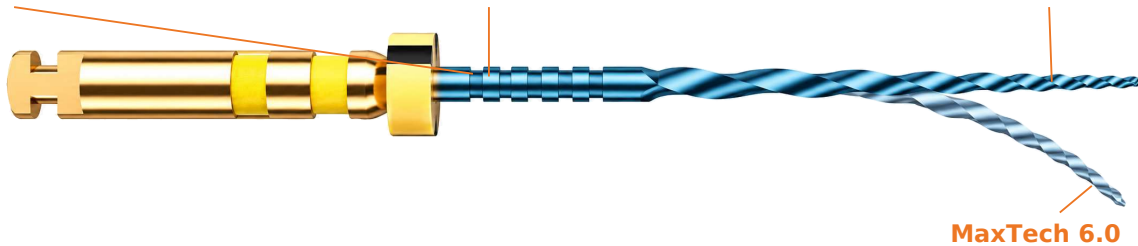

File sequence guide

SCOUTING

U-BLUE

WHEN NEEDED



File	SS K	SS K	SX	S1	S2	F1	F2	F3	F4	F5
Size / taper	#10	#15	019/.04	017/.02	020/.04	020/.07	025/.08	030/.09	040/.06	050/.05
Speed (rpm)	-	-	400	400	400	400	400	400	400	400
Torque (N·cm)	-	-	3.0	3.0	3.0	3.0	3.0	3.0	3.0	3.0

K #10 → K #15 → SX → S1 → S2 → F1 → F2 → F3. F4 / F5 when needed.

Settings transcribed from the supplied NIC guide. Follow the manufacturer's instructions for use.

CANAL DRYING



[Click to Visit Website](#)



Absorbent paper points

Accurate, absorbent and length marked for convenient reference.

ISO paper points

For U-Series Files



ISO 15-40 | 200 points



F1-F3 | 100 points

Designed for dependable drying

High solidity & absorbency

Reliable performance and ease of use.

Accurate

Precision for routine canal drying.

Colour-coded

Easy recognition when selecting points.

Consistent quality

Uniform quality across the range.

Clean, lintless surface

Gamma-sterilized to support infection control.

Convenient one-hand sliding-open package.

Both ISO and U-Series paper points are length marked.

OBTURATION



[Click to Visit Website](#)



Gutta-percha points

Machine-made precision. Both options are length marked.

ISO gutta-percha

For U-Series Files



ISO 15-40 | 120 points



VTF F1-F3 | 60 points

Individual point: For illustration purposes only.

01 Machine-made

Uniform taper and shape, with improved apical filling capacity compared with handmade points. Designed to match the corresponding reamer or file.

02 Precise

Strict manufacturing standards for quality and a precise fit.

03 Consistent quality

Designed to match both constant- and variable-taper preparations.

RADIOPAQUE / ELASTIC / NON-LATEX

Visible on X-rays. Flexible during insertion and compaction.

Depth markings: 18, 19, 20 and 22 mm; illustrated positions are approximate.

ORDER THROUGH UNIVERSAL DENTAL



[Click to Visit Website](#)



Ordering information

Select your product, size, file length and pack quantity.

PRODUCT	LENGTH	AVAILABLE SIZES / OPTIONS	PACK
H-Files Stainless Steel	21 mm 25 mm	8, 10, 15, 20, 25, 30, 35, 40 Assorted: 15-40 or 45-80	6 pcs
K-Files Stainless Steel	21 mm 25 mm	6, 8, 10, 15, 20, 25, 30, 35, 40 Assorted: 15-40 or 45-80	6 pcs
U-Blue files NiTi rotary files	SX: 19 mm Others: 21 / 25 mm	SX, S1, S2, F1, F2, F3 Assorted: SX-F3	6 pcs
ISO paper points 2% taper · Length marked	-	15, 20, 25, 30, 35, 40 Assorted: 15-40 or 45-80	200 pcs
U-Series Paper Points Tapered Length Marked	-	F1, F2, F3	100 pcs
ISO gutta-percha 2% taper · Length marked	-	15, 20, 25, 30, 35, 40 Assorted: 15-40 or 45-80	120 pcs
U-Series Gutta-percha Tapered Length Marked	-	F1, F2, F3	60 pcs

ORDER NOW

WHATSAPP

0300-0494959

Universal Dental

WEBSITE

www.universaldental.com.pk

FDA/CE & DRAP Registered